

Anabolic Lab™

Analytical Lab Testing for
Anabolic Steroid Harm Reduction

<https://anaboliclab.com>

Lab Report – April 27, 2016

Geneza Pharmaceuticals
GP Oxy (50mg/tablet)
Oxymetholone



If you like the attached report, please visit [our website](https://anaboliclab.com) to find out how you can support independent and objective lab testing.



Simec AG

Areal Bleiche West Postfach 413 CH - 4800 Zofingen
Telefon +41 62 752 83 08 Fax +41 62 752 83 09
E-Mail info@simec.ch www.simec.ch

Test report

Number **C 42759**

Object of analysis: **Geneza GP Oxy**

Customer: **AnabolicLab.com**

Samples: **1**

Batch: **GP672H, Exp. 09/2017**

Receipt: **01.04.2016**

Subject: **Assay, Identification Oxymetholone**

Method(s): **HPLC**

Sample / Analysis:	Method	Status	Result	Specification	complies
<u>Oxymetholone:</u>					
Identification	HPLC	N	complies	complies	yes
Content	HPLC	N	50.9 mg/tablet	50 mg/tablet	-

Method-status: G = GMP A = accredited V = generally validated P = validated on product N = not validated
E = external lab

Remarks: The sample(s) mentioned above have been analysed as they have been sent to us by the client. We have not controlled external processes like manufacturing, labelling, sampling, shipping, storage. These results are for information only and do not compensate for correct quality control by the manufacturer/distributor. Generally, pharmaceutical products have to be produced and distributed under full GMP/GDP regime, including GMP-compliant analyses with validated methods.

Signatures: **controlled & released:**

Staff:

completed:

27.04.2016


The enclosed results refer exclusively to the object of examination described above. The measuring accuracy is available on request. Only handwritten signatures are valid. Without written permission, it is not allowed to publish single parts of this report.

page 1 of 1

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 24.04.2016 19:41:21
 Vial: 42
 Probenidentifikation: C42759 A
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 85, 250 x 4,6, C18
 Detektion: UV 240nm

Druckdatum: 25.04.2016 09:00:03
 Injektvolumen: 20

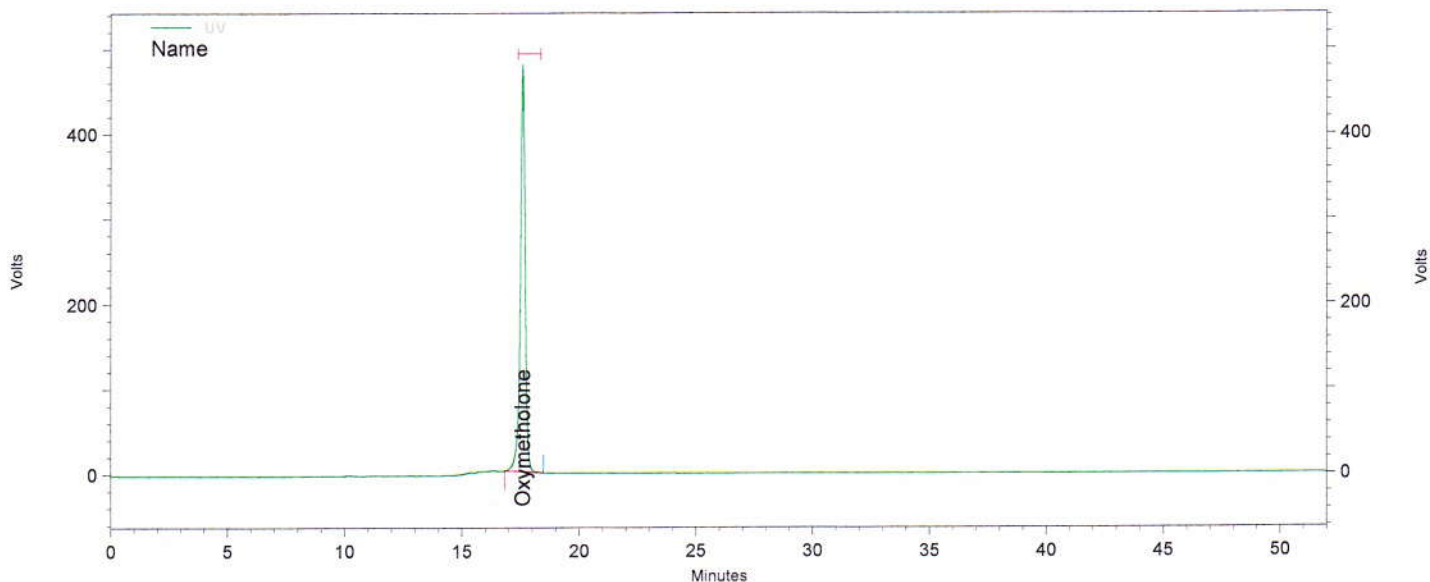

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone	17,597	27131388	0,000	100,000

Totals		27131388	0,000	100,000
--------	--	----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\April 2016\Diverse
 Steroide\19,04,16\054_C42759 A-Rep1.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_UV290.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 24.04.2016 20:34:43 Druckdatum: 25.04.2016 09:00:07
 Vial: 42 Injektvolumen: 20
 Probenidentifikation: C42759 A
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 85, 250 x 4,6, C18
 Detektion: UV 240nm

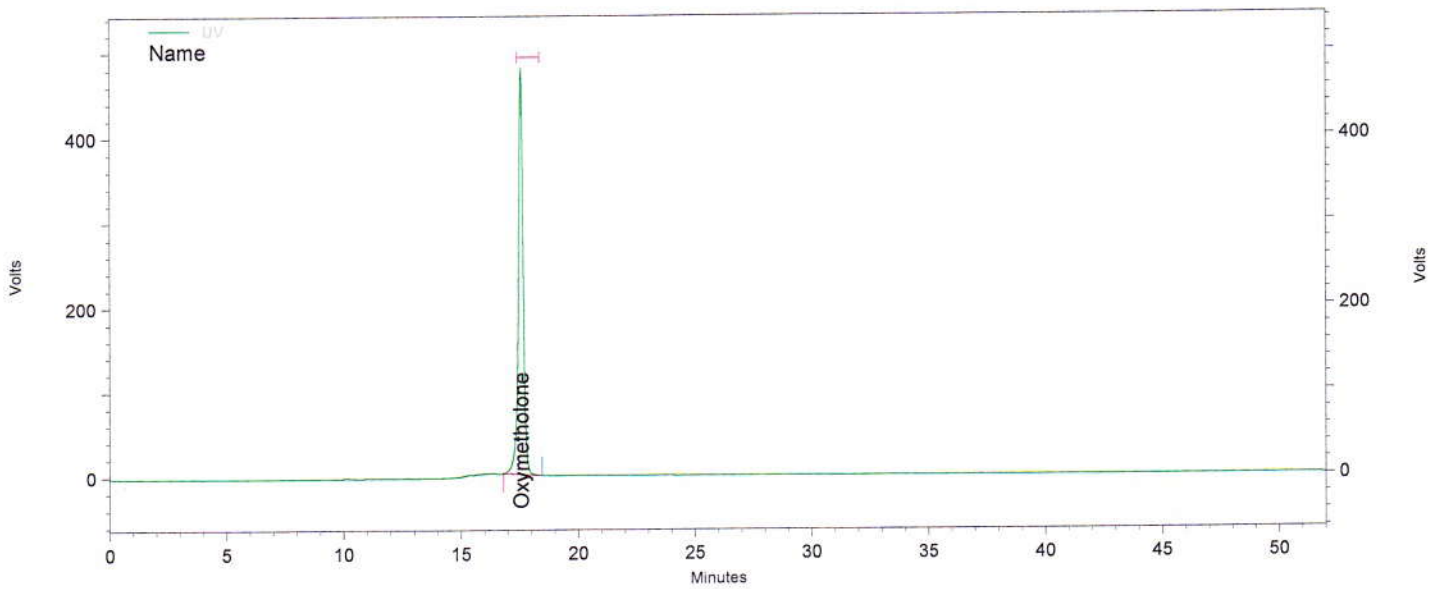

Channel A


UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	17,563	27172004	0,000	100,000
Totals		27172004	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\April 2016\Diverse
 Steroide\19,04,16\054_C42759 A-Rep2.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_UV290.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 24.04.2016 21:28:03 Druckdatum: 25.04.2016 09:00:11
 Vial: 42 Injektvolumen: 20
 Probenidentifikation: C42759 A
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 85, 250 x 4,6, C18
 Detektion: UV 240nm

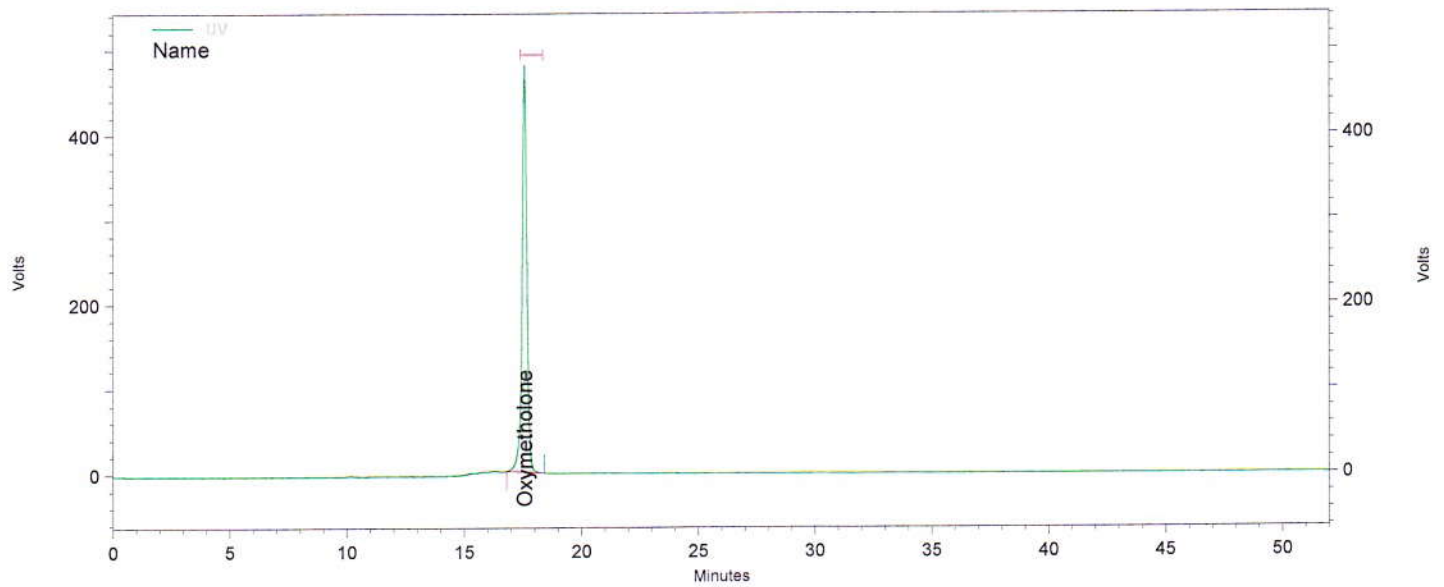

Channel A


UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	17,577	27150790	0,000	100,000
Totals		27150790	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\April 2016\Diverse Steroide\19,04,16\054_C42759 A-Rep3.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_UV290.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 24.04.2016 22:21:25 Druckdatum: 25.04.2016 09:00:15
 Vial: 43 Injektvolumen: 20
 Probenidentifikation: C42759 B
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 85, 250 x 4,6, C18
 Detektion: UV 240nm

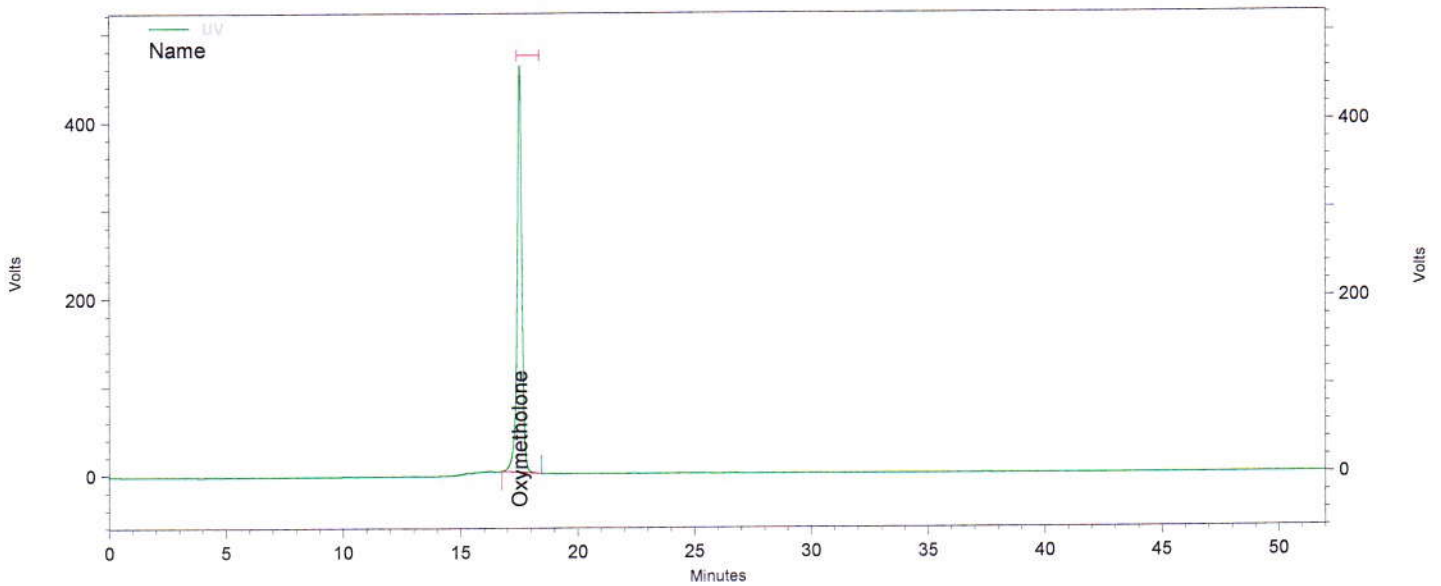

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone	17,523	26115025	0,000	100,000
Totals		26115025	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\April 2016\Diverse
 Steroide\19,04,16\055_C42759 B-Rep1.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_UV290.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 24.04.2016 23:14:52
 Vial: 43
 Probenidentifikation: C42759 B
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 85, 250 x 4,6, C18
 Detektion: UV 240nm

Druckdatum: 25.04.2016 09:00:19
 Injektvolumen: 20

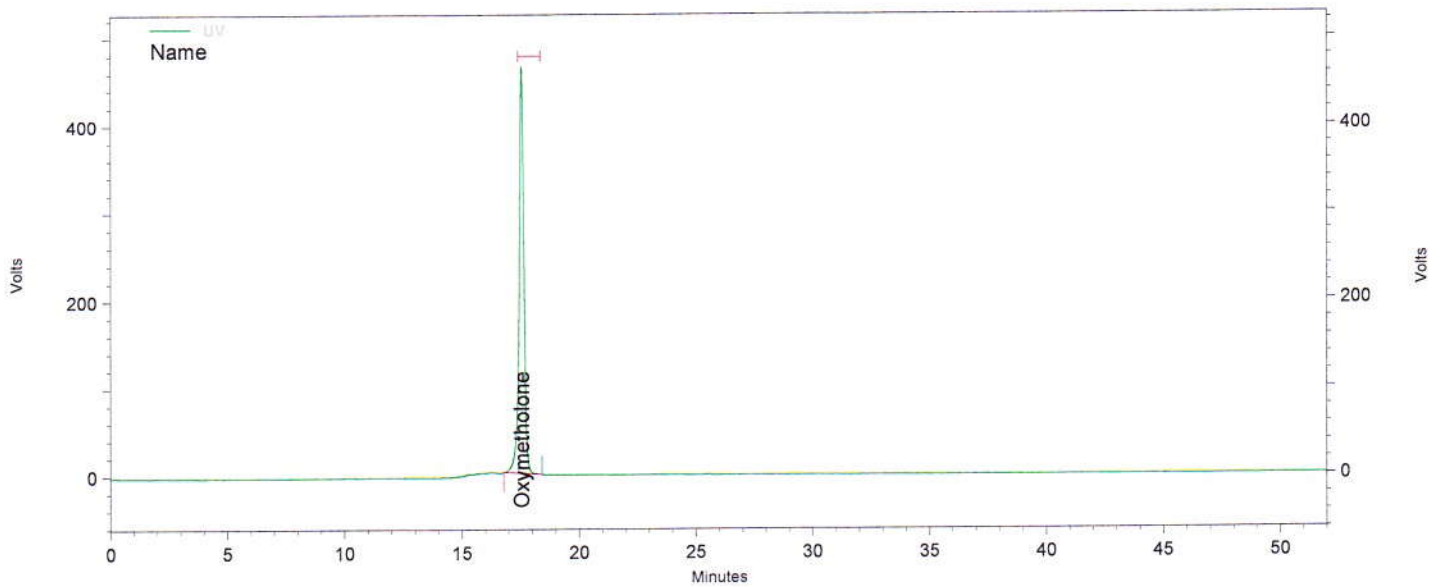

Channel A

UV Results


Name	Time	Area	Amount	Area Percent
Oxymetholone	17,553	26319373	0,000	100,000
Totals		26319373	0,000	100,000

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\April 2016\Diverse
 Steroide\19,04,16\055_C42759 B-Rep2.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_UV290.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 25.04.2016 00:08:21
 Vial: 43
 Probenidentifikation: C42759 B
 Operator: admin
 Chrom.-Beding.: ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule: Nr. 85, 250 x 4,6, C18
 Detektion: UV 240nm

Druckdatum: 25.04.2016 09:00:23
 Injektvolumen: 20

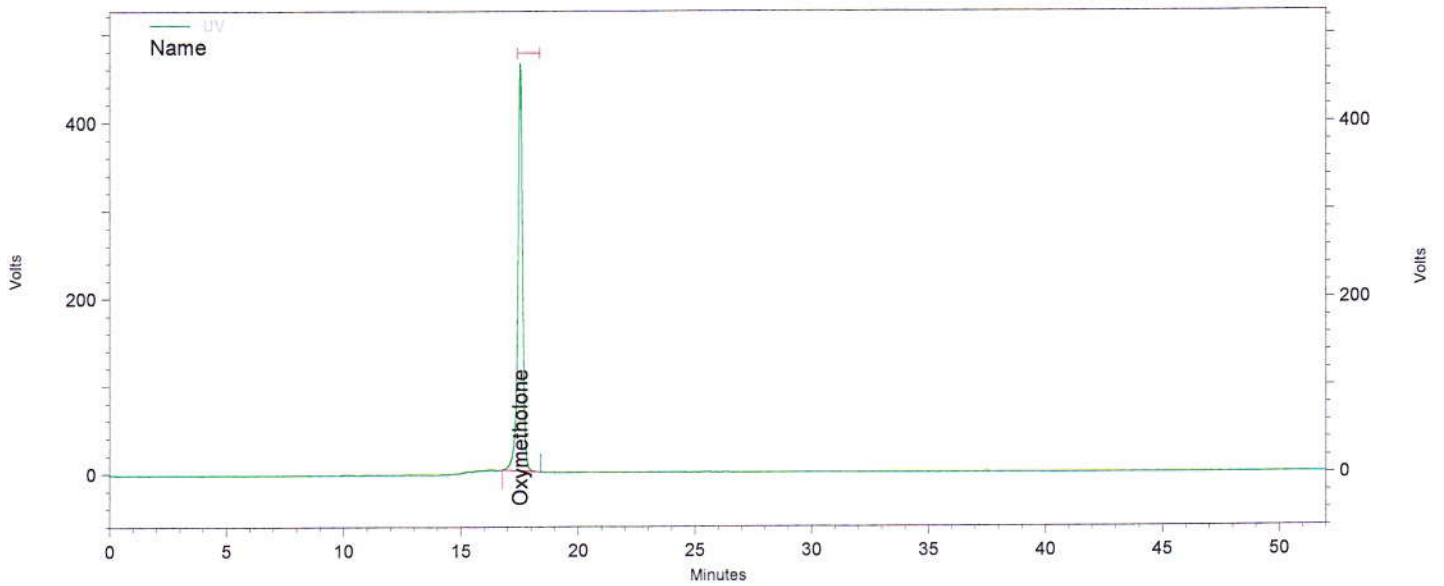

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone	17,530	26159936	0,000	100,000


Totals		26159936	0,000	100,000
--------	--	----------	-------	---------

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\April 2016\Diverse
 Steroide\19,04,16\055_C42759 B-Rep3.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_UV290.met

Report steroid analysis

HPLC: Nr.5
 Analysendatum: 25.04.2016 01:01:50
 Vial: 1
 Probenidentifikation: Blank
 Operator: admin
 Chrom.-Beding. ACN / Wasser, gradient
 Temperatur: 40 °C
 Säule Nr. 85, 250 x 4,6, C18
 Detektion: UV 240nm

Druckdatum: 25.04.2016 09:00:27
 Injektvolumen: 20

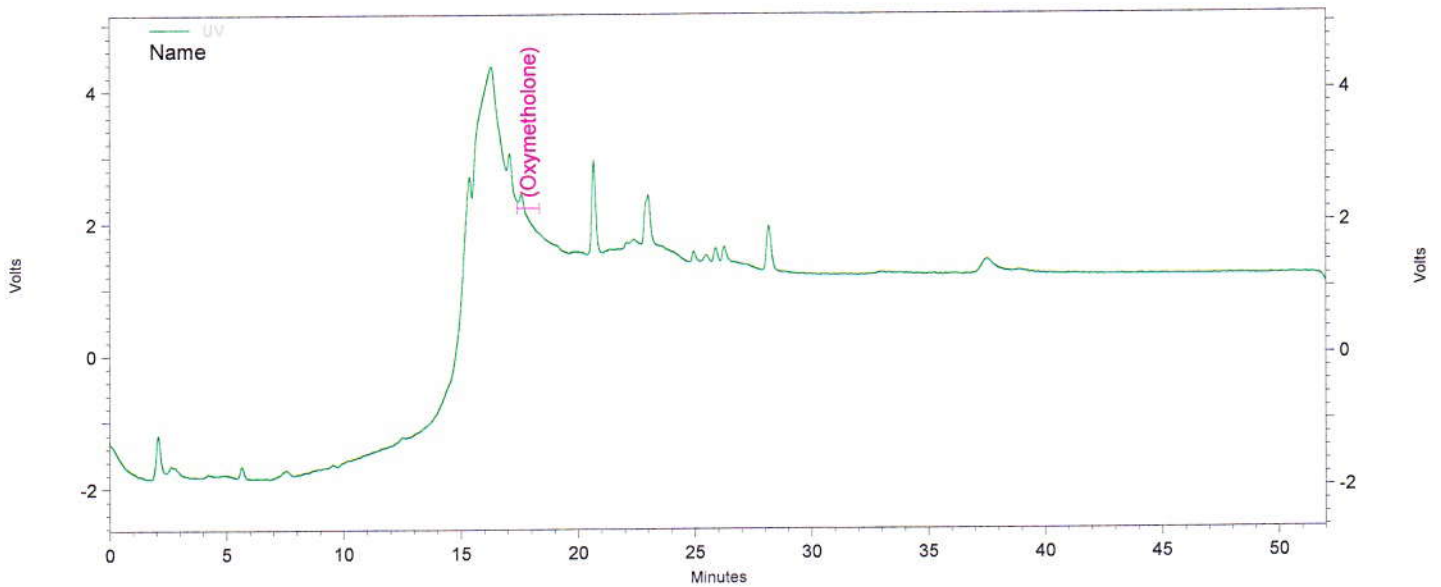

Channel A

UV Results

Name	Time	Area	Amount	Area Percent
Oxymetholone			0,000 BDL	
Totals				

Daten File: C:\EZChrom Elite\Enterprise\Projects\Simec\Data\2016\April 2016\Diverse Steroide\19,04,16\056_Blank.dat
 Methode: C:\EZChrom Elite\Enterprise\Projects\Simec\Method\Steroide\Oxymetholone_UV290.met